

Product brochure

05.09.2026

KERN PLJ 1200-3A

Range of precision balances with enormous weighing capacities – ideal for heavy tare containers or large samples


ENGLISH



CONTACT

KERN & SOHN GmbH
Ziegelei 1
72336 Balingen
Germany

+49 7433 9933-0
info@kern-sohn.com
www.kern-sohn.com

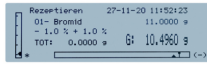

FEATURES

- Convenient recipe-weighing with the recipe database, up to 99 recipes can be stored, each with up to 20 recipe ingredients with name and target value
- Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display
- Rapid and efficient operation thanks to the graphics display
- Dosage aid: High stability mode and other filter settings can be selected
- Simple, clear user interface on the display in the following languages: DE, EN, FR, IT, ES, PT
- Automatic internal adjustment, guarantees high degree of accuracy and makes the balance independent of its location of use
- Ring-shaped draught shield standard, weighing space $\varnothing \times H$ 150×60 mm
- Protective working cover included with delivery

TECHNICAL DATA

- Backlit LCD graphic display, digit height 15 mm
- Dimensions weighing surface $\varnothing$ 110 mm, stainless steel
- Overall dimensions
W×D×H 210×340×120 mm (without draught shield)
W×D×H 210×340×160 mm (with draught shield)
- Weighing principle: Force compensation
- Permissible ambient temperature 15 °C/35 °C

MAIN SCOPE OF APPLICATIONS

- Laboratory applications
- Research and development
- Jewellery and precious metal testing
- Environmental analysis
- Electronics industry

PICTOGRAMS

Standard



Option



MODEL OVERVIEW

Model	Weighing capacity [Max]	Readability [d]	Reproducibility	Linearity	Dimensions weighing surface
PLJ 420-3F	420 g	0,001 g	0,001 g	± 0,003 g	110 mm
PLJ 720-3A	720 g	0,001 g	0,001 g	± 0,002 g	110 mm
PLJ 1200-3A	1200 g	0,001 g	0,001 g	± 0,003 g	110 mm
PLJ 2000-3A	2100 g	0,001 g	0,002 g	± 0,004 g	110 mm
PLJ 4200-2F	4200 g	0,01 g	0,02 g	± 0,04 g	160 mm
PLJ 6200-2A	6200 g	0,01 g	0,01 g	± 0,05 g	160 mm

ENGLISH

ACCESSORIES

Model	Description
315-111-206	Safety set 315-111-206, included weight 50 g (F1); 1000 g (E2)
572-926	RS232 interface cable 572-926
ALT-A02	Density determination set ALT-A02
PLJ-A01S05	Protective working cover PLJ-A01S05, scope of delivery: 5 items
PLJ-A02	Hook for underfloor weighing PLJ-A02
PLJ-A09	Protective working cover PLJ-A09
PLJ-A09S05	Protective working cover PLJ-A09S05, scope of delivery: 5 items
YKC-01	Thermal printer YKC-01
YKG-01	Dot Matrix Printer YKG-01
YKH-01	Thermal printer YKH-01
YKI-02	RS232-Bluetooth-Adapter YKI-02
YKI-10	RS232-Ethernet adapter YKI-10 (RS-232 standard, Ethernet, RS-485)
YKI-11	RS232-Ethernet adapter YKI-11 (Ethernet, RS-232 standard - 2x)
YKI-12	RS232-WiFi adapter YKI-12
YKT-01	Needle printer YKT-01
YPS-03	Antivibration table YPS-03

Model	Description
YPS-04	Anti-vibration plate KERN YPS-04 (400×450×60 mm, 30 kg)
YPS-05	Anti-vibration plate KERN YPS-05 (565×450×60 mm, 42 kg)

SOFTWARE

Model	Description
SCD-4.0-DL	Software SCD-4.0-DL (Scope of supplies: Download link for 1 license)
SCD-4.0-DLS05	Software SCD-4.0-DLS05 (Scope of supplies: Download link for 5 licenses)
SCD-4.0-PRO-DL	Software SCD-4.0-PRO-DL (Scope of supplies: Download link for 1 license)
SCD-4.0-PRO-DL-AS05	Software SCD-4.0-PRO-DL-AS05 (Scope of supplies: Download link for 5 licenses)

SERVICES

Model	Description
963-103	DAkkS calibration certificate 963-103
969-517	Confirmation of the manufacturer's specification (only in combination with DAkkS calibration certificate) 969-517
970-107	Warranty extension (+2 years) 970-107

KERN Pictogram Overview

	Internal adjustment: Adjusting accuracy using internal adjustment weight
	Data interface RS-232: For connection to printers, PCs or other peripheral devices for data transfer over longer distances. Network in bus topology possible
	GLP/ISO protocol: The measuring device generates a GLP-compliant printout
	Value & time: The measuring device displays the value, date, and time, regardless of the connected printer
	Piece counting: Reference quantities selectable. Display can be switched from piece to weight
	Formulation: Internal memory for complete recipes with name and target value of recipe components. Display-supported user guidance
	Density determination: The density of liquids and solids with a density of $\leq / \geq 1$ is determined directly in the measuring device
	Percentage determination: determining the deviation in % from the target value (100 %)
	Measuring with tolerance range (checkweighing): Upper and lower limit values can be programmed, e.g. for sorting, portioning or for measuring coatings or material thicknesses. The process is supported by an acoustic and/or optical signal
	Suspended weighing: Load support with hook on the underside of the balance
	Universal plug-in power supply: With universal input and optional input socket adapters for EU, CH, GB or EU, CH, GB, US or EU, CH, GB, US, AUS
	Weighing principle electromagnetic force compensation: Coil inside a permanent magnet. For the most accurate weighings
	Package shipment by courier service: The duration of internal product provision is one day
	Accredited calibration (DKD): The duration of the accredited calibration is 3 working days